



ADT Interface Specifications

Product Name: HistoTrac

Component Name: SystemLink ADT Service

Version 2.0

CONFIDENTIAL TO SYSTEMLINK, INC. AND ITS CLIENTS

Table of Contents

1	INTRODUCTION	1
1.1	BUSINESS CASE	1
1.2	SOFTWARE CHANGE INTENDED USE	1
2	GENERAL WORKFLOW	2
3	HISTOTRAC AND 3RD PARTY SYSTEM COMMUNICATION	3
3.1	COMMUNICATION PROTOCOL AND TRANSACTION FRAMING	3
3.2	ACKNOWLEDGMENTS (ACK/NAK)	3
4	HISTOTRAC HL7 ADT TRANSACTION AND SEGMENTS	4
4.1	SUPPORTED ADT EVENTS	4
4.2	ADT MESSAGE	4
4.2.1	MSH – Message Header	4
4.2.2	EVN – Event Segment	5
4.2.3	PID – Patient Identification	5
4.2.4	PV1 – Patient Visit	6
4.2.5	NTE – Notes and Comments	7
4.3	MESSAGE ACKNOWLEDGEMENT	7
4.3.1	MSH – Message Header	7
4.3.2	MSA – Message Acknowledgement	8

1 Introduction

This specification details the structure for an inbounds ADT interface. This interface allows a third party system to send an Admission, Discharge, and Transfer (ADT) data to HistoTrac that will create or update patient data in HistoTrac. The ADT data is recorded in HistoTrac and are placed in a "holding" database until the Medical Record Number is used in HistoTrac. If the medical record number is referenced in HistoTrac the corresponding ADT data is copied into the HistoTrac repository. Unused ADT data is purged after a specified number of days.

1.1 Business Case

Third party systems that generate ADT data can pass this information into the HistoTrac application. HistoTrac users are currently manually entered into HistoTrac. Since this data has already been placed in a third party system it is more efficient to pass the ADT data to HistoTrac through an automated interface. By utilizing the interface the ADT data will no longer need to be entered into HistoTrac reducing transcription errors and delays in getting paperwork to the laboratory. ADT data received through the interface in HistoTrac will be placed in a temporary database until the Medical Record number is used in the application. Data that is not referenced in HistoTrac within a specified number of days is removed to reduce the volume of data in the ADT database.


1.2 Software Change Intended Use

The ADT Interface is developed so that third party software systems can send patient information directly to the HistoTrac repository.

2 General Workflow

With the ADT interface the third party system will send a correctly formatted HL7 ADT transaction to a specific IP Address and port. The HistoTrac HL7 Manager will be listening for ADT traffic on the specified IP address and port. When a message is received an acknowledge message is created and is transmitted back to the sending IP Address and port. This acknowledgement transaction signals the remote system that HistoTrac is ready to accept a new ADT transaction.

Data Flow

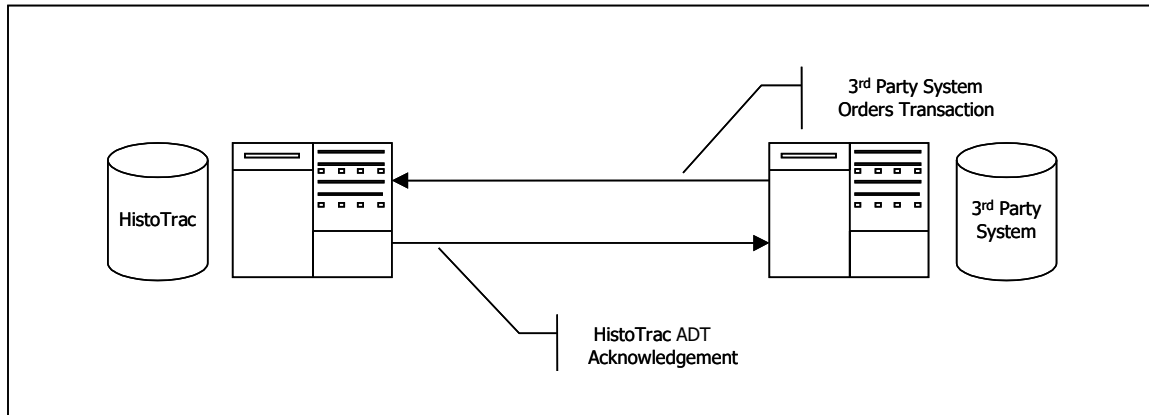

3 HistoTrac and 3rd Party System Communication

3.1 Communication Protocol and Transaction Framing

All transactions will be sent using TCP/IP as the communication protocol. With the reliability of TCP, no block check calculations or data length counts are necessary. The following applies to the communication envelope for all messages. The characters listed in <> below represent ASCII characters.

- All messages and acknowledgments will begin with <11>
- All segments will end with <13>
- All messages and acknowledgments will terminate with <28><13>

3.2 Acknowledgments (ACK/NAK)

The receiving system (HistoTrac/server) must return a communication level acknowledgment to the sending system (3rd party system/client) for every message that is received. Once the client has received the acknowledgment, it may then send the next message. For example, when an ADT transaction is sent to HistoTrac, HistoTrac will return an acknowledgment to the originating system for that message on the same logical socket where it was received. The other system can then send the next ADT transaction to HistoTrac.

Acknowledgments can be either positive (ACK) or negative (NAK). For all ACK/NAK messages, the HL7 segments MSH and MSA must be returned. The MSA; 1 is valued with "CA" for ACK and with "CE" for NAK.

4 HistoTrac HL7 ADT Transaction and Segments

4.1 Supported ADT Events

HL7 (v2.3) Triggers	Functionality	Received By HistoTrac
A01	Admit Inpatient	Y
A02	Transfer in or outpatient	Y
A03	Discharge in or outpatient	Y
A04	Register outpatient or outside patient	Y
A05	Preadmit patient	N
A06	Transfer outpatient to inpatient	Y
A07	Transfer inpatient to outpatient	Y
A08	Demographic update	Y
A11	Cancel admit	N
A13	Cancel discharge	Y
A17	Swap patients	N
A23	Delete a patient record	N
A28	Add person information	N
A29	Delete person information	N
A31	Update person information	N
A34	Merge patient information	Y
A35	Merge account information	N

4.2 ADT Message

The ADT transaction is separated into the following five segments:

MSH – Message Header
 EVN – Event Information
 PID – Patient Identification
 PV1 – Patient Visit
 [{NTE}] – Notes and Comments

The following table outlines the contents of each segment. Each segment is required in the ADT transaction.

Legend
R – Required
O – Optional
NS – Not Supported

4.2.1 MSH – Message Header

Sequence #	Name	Max Length	Field Requirement	Notes
1	Field Separator	1	R	
2	Encoding Characters	4	R	component delimiter (^) repeat delimiter (~) escape character (\) subcomponent delimiter (&)
3	Sending Application	50	R	
4	Sending Facility	50	O	
5	Receiving Applications	50	O	
6	Receiving Facility	50	O	
7	Date/Time of Message	12	R	CCYYMMDDHHMM
8	Security	N/A	NS	
9	Message Type	7	R	Contains message type and trigger event.

Information contained in this document is confidential to SystemLink, Inc and its principals and may not be used without permission.

10	Message Control ID	15	R	This sequence is valued with a number that is unique to that message. The receiving system will echo this value in the MSA message.
11	Processing ID	3	R	T for Training or P for Production
12	Version ID	8	R	Specifies the HL7 version in use
13	Sequence Number	N/A	NS	
14	Continuation Pointer	N/A	NS	
15	Accept Acknowledgement Type	N/A	NS	
16	Application Acknowledgement Type	N/A	NS	
17	Country Code	N/A	NS	
18	Character Set	N/A	NS	
19	Principal Language of Message	N/A	NS	

4.2.2 EVN – Event Segment

Sequence #	Name	Max Length	Field Requirement	Notes
1	Event Type Code	5	R	
2	Date/Time of Event	12	R	CCYYMMDDHHMM
3	Date/Time Planned Event	N/A	NS	
4	Event Reason Code	N/A	NS	
5	Operator ID	N/A	NS	

4.2.3 PID – Patient Identification

Sequence #	Name	Max Length	Field Requirement	Notes
1	Set ID – Patient ID	N/A	NS	
2	Patient ID (External ID)	N/A	NS	
3	Patient ID (Internal ID)	20	R	Component 1 = Medical Record Number Component 2 = "" Component 3 = MulHos Code Example: Z442345^^FZQ
4	Alternate Patient ID	20	O	
5	Patient's Name	55	R	Component 1=LAST NAME Component 2=FIRST NAME Component 3=MIDDLE INITIAL Example: SMITH^JOHN^D
6	Mother's Maiden Name	40	O	
7	Date of Birth	8	O	CCYYMMDD
8	Sex	15	O	
9	Patient Alias	N/A	NS	
10	Race	15	O	
11	Patient Address	40, 40, 30, 2, 10	O	Component 1=address line 1 Component 2=address line 2 Component 3=city Component 4=state Component 5=zip code Example: 123 Oak Street ^ Suite 100 ^ Anywhere ^ VA ^ 12345
12	Country Code	15	O	
13	Phone Number – Home	15	O	
14	Phone Number –	15	O	

Information contained in this document is confidential to SystemLink, Inc and its principals and may not be used without permission.

	Business			
15	Language – Patient	N/A	NS	
16	Marital Status	N/A	NS	
17	Religion	N/A	NS	
18	Patient Account Number	25	O	
19	SSN Number – Patient	15	O	
20	Driver's License – Patient	N/A	NS	
21	Mother's Identifier	N/A	NS	
22	Ethnic Group	15	O	
23	Birth Place	N/A	NS	
24	Multiple Birth Indicator	N/A	NS	
25	Birth Order	N/A	NS	
26	Citizenship	N/A	NS	
27	Veteran's Military Status	N/A	NS	
28	Nationality	N/A	NS	
29	Patient Death Date/Time	N/A	NS	
30	Patient Death Indicator	N/A	NS	

4.2.4 PV1 – Patient Visit

Sequence #	Name	Max Length	Field Requirement	Notes
1	Set ID PV1	N/A	NS	
2	Patient Class	1	R	
3	Assigned Patient Location	20, 20, 20	O	Component 1 = Unit Component 2 = Room Component 3 = Bed Example: unit^room^bed
4	Admission Type	N/A	NS	
5	Pre-Admit Number	N/A	NS	
6	Prior Patient Location	N/A	NS	
7	Attending Doctor	10, 25, 25	O	Component 1= Physician Code Component 2 = Last Name Component 3 = First Name Example: GRT1^GRANT^JOHN
8	Referring Doctor	N/A	NS	
9	Consulting Doctor	N/A	NS	
10	Hospital Service	N/A	NS	
11	Temporary Location	N/A	NS	
12	Pre-Admit Test Indicator	N/A	NS	
13	Re-Admission Indicator	N/A	NS	
14	Admit Source	N/A	NS	
15	Ambulatory Status	N/A	NS	
16	VIP Indicators	N/A	NS	
17	Admitting Doctor	10, 25, 25	O	Component 1= Physician Code Component 2 = Last Name Component 3 = First Name Example: GRT1^GRANT^JOHN
18	Patient Type	N/A	NS	
19	Visit Number	15, 15	R	Component 1 = Episode ID Component 2 = Event ID Example: 1^2
20	Financial Class	N/A	NS	
21	Charge Price Indicator	N/A	NS	
22	Courtesy Code	N/A	NS	
23	Credit Rating	N/A	NS	

Information contained in this document is confidential to SystemLink, Inc and its principals and may not be used without permission.

24	Contract Code	N/A	NS	
25	Contract Effective Date	N/A	NS	
26	Contract Amount	N/A	NS	
27	Contract Period	N/A	NS	
28	Interest Code	N/A	NS	
29	Transfer to Bad Debt Code	N/A	NS	
30	Transfer to Bad Debt Date	N/A	NS	
31	Bad Debt Agency Code	N/A	NS	
32	Bad Debt Transfer Amount	N/A	NS	
33	Bad Debt Recovery Amount	N/A	NS	
34	Delete Account Indicator	N/A	NS	
35	Delete Account Date	N/A	NS	
36	Discharge Disposition	N/A	NS	
37	Discharged To Location	N/A	NS	
38	Diet Type	N/A	NS	
39	Servicing Facility	N/A	NS	
40	Bed Status	N/A	NS	
41	Account Status	N/A	NS	
42	Pending Location	N/A	NS	
43	Prior Temporary Location	N/A	NS	
44	Admit Date/Time	8	R	CCYYMMDD
45	Discharge Date/Time	8	O	CCYYMMDD
46	Current Patient Balance	N/A	NS	
47	Total Charges	N/A	NS	
48	Total Adjustments	N/A	NS	
49	Total Payments	N/A	NS	
50	Alternate Visit ID	N/A	NS	
51	Visit Indicator	N/A	NS	
52	Other Health Care Provider	N/A	NS	

4.2.5 NTE – Notes and Comments

Sequence #	Name	Max Length	Field Requirement	Notes
1	Set ID – NTE	4	O	
2	Source Of Comment	8	O	
3	Comment	2000	O	

4.3 Message Acknowledgement

The Message Acknowledgement transaction is separated into the following two segments:

- 1 MSH – Message Header
- 2 MSA – Message Acknowledgement

The following table outlines the contents of each segment. Each segment is required in the ADT transaction.

Legend
R – Required
O – Optional
NS – Not Supported

4.3.1 MSH – Message Header

Sequence #	Name	Max Length	Field Requirement	Notes
1	Field Separator	1	R	
2	Encoding Characters	4	R	component delimiter (^) repeat delimiter (~) escape character (\) subcomponent delimiter (&)
3	Sending Application	50	R	
4	Sending Facility	50	O	
5	Receiving Applications	50	O	
6	Receiving Facility	50	O	
7	Date/Time of Message	12	R	CCYYMMDDHHMM
8	Security	N/A	NS	
9	Message Type	7	R	Contains message type and trigger event.
10	Message Control ID	15	R	This sequence is valued with a number that is unique to that message. The receiving system will echo this value in the MSA message.
11	Processing ID	3	R	T for Training or P for Production
12	Version ID	8	R	Specifies the HL7 version in use
13	Sequence Number	N/A	NS	
14	Continuation Pointer	N/A	NS	
15	Accept Acknowledgement Type	N/A	NS	
16	Application Acknowledgement Type	N/A	NS	
17	Country Code	N/A	NS	
18	Character Set	N/A	NS	
19	Principal Language of Message	N/A	NS	

4.3.2 MSA – Message Acknowledgement

Sequence #	Name	Max Length	Field Requirement	Notes
1	Acknowledgement Code	2	R	Indicates whether positive (ACK) or negative (NAK) acknowledgment. CA = ACK CR = NAK
2	Message Control ID	15	R	The number sent in MSH; 10 of the inbound message is echoed back to the client in this sequence.
3	Text Message	100	O	This field may be valued but only if returning a NAK and only under certain error conditions. A description of the error may be included here for troubleshooting purposes
4	Expected Sequence Number	N/A	NS	HistoTrac does not use sequence number protocol
5	Delayed ACK Type	N/A	NS	
6	Error Condition	N/A	NS	